
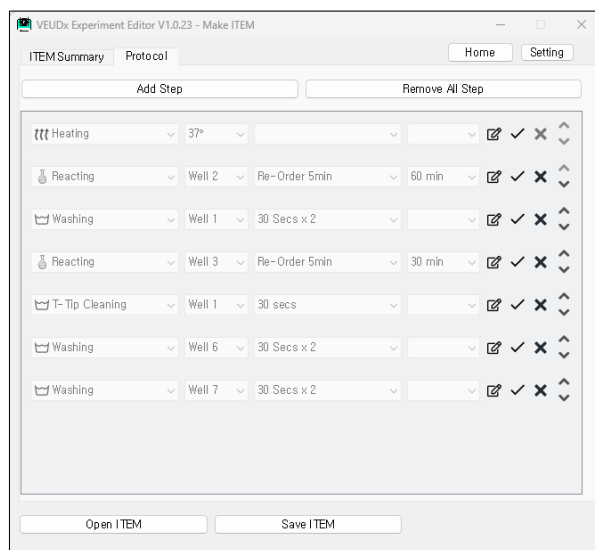


# Operation Manual

## VEUDx Experiment Editor

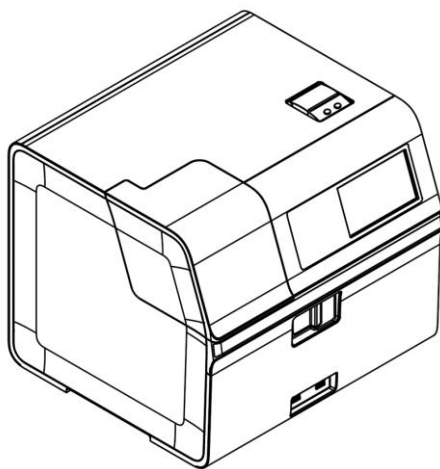

UM-VEUDx-1.0

## Index

|                                      |    |
|--------------------------------------|----|
| 1. Installation.....                 | 3  |
| 2. Term definition.....              | 5  |
| 2.1 ITEM .....                       | 5  |
| 2.2 LOT .....                        | 5  |
| 2.3 QC Material LOT.....             | 5  |
| 3. Start.....                        | 6  |
| 3.1 Start screen.....                | 6  |
| 4. Make ITEM .....                   | 6  |
| 4.1 ITEM creation start screen.....  | 6  |
| 4.2 Edit ITEM Summary .....          | 8  |
| 4.3 Edit Protocol .....              | 10 |
| 4.4 Save ITEM and Installation ..... | 11 |
| 5. Make LOT .....                    | 13 |
| 5.1 Select ITEM file .....           | 13 |
| 5.2 Make LOT .....                   | 13 |
| 6. Make QC Material LOT.....         | 17 |
| 6.1 Select ITEM file .....           | 17 |
| 6.2 Make QC Material LOT.....        | 17 |
| 6. Settings.....                     | 19 |
| 6.1 Open Settings.....               | 19 |
| 6.2 Unit.....                        | 19 |
| 6.3 Development menu.....            | 19 |

# 1. Installation


## 1.1 Installation

### 1.1.1 Execution

Execute VEUDx Experiment Editor Setup V1.x.x.exe on PC.


### 1.1.2 Language selection

Choose the language you want to use

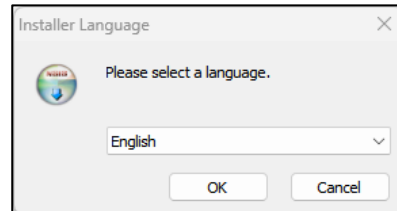


### 1.1.3 Start installation

Check the contents and click 'Next'.

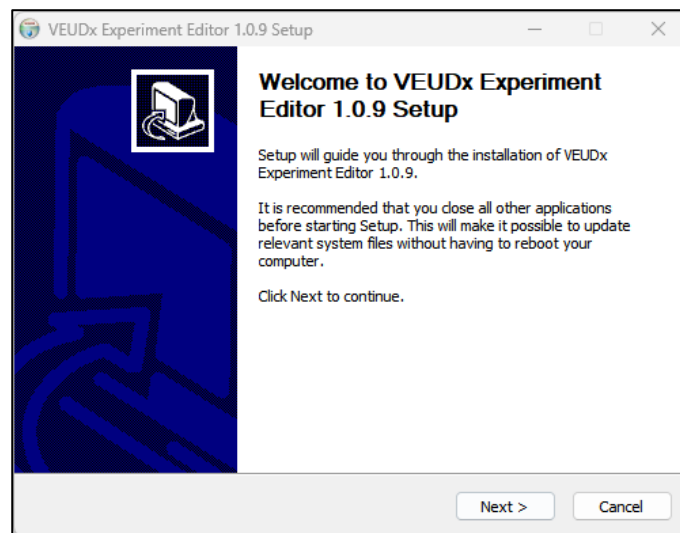


### 1.1.4 Check License

Read the license agreement and click 'I Agree'.

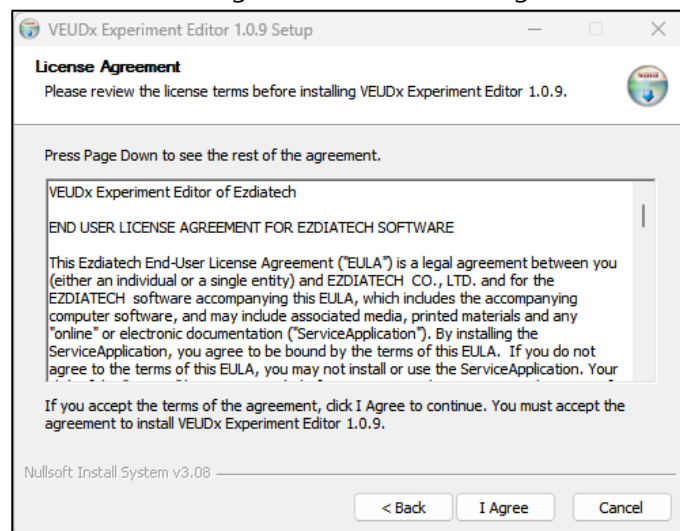


### 1.1.5 Installation location

Check the installation path and click 'Install'.

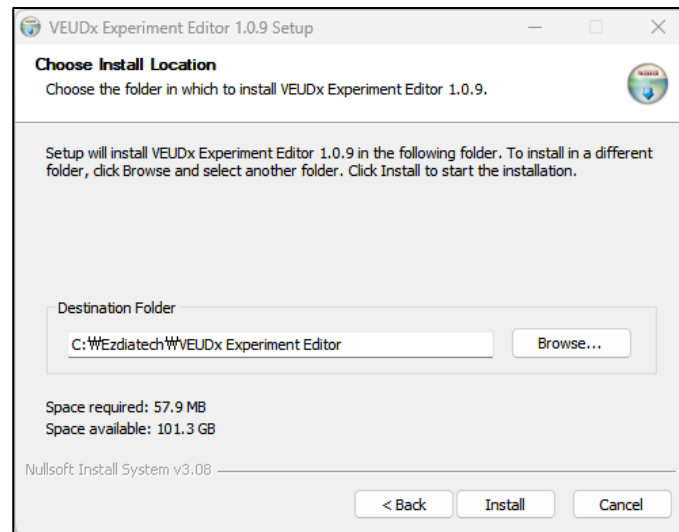

### 1.1.6 Installation completion

Check to run the program and click 'Finish'.

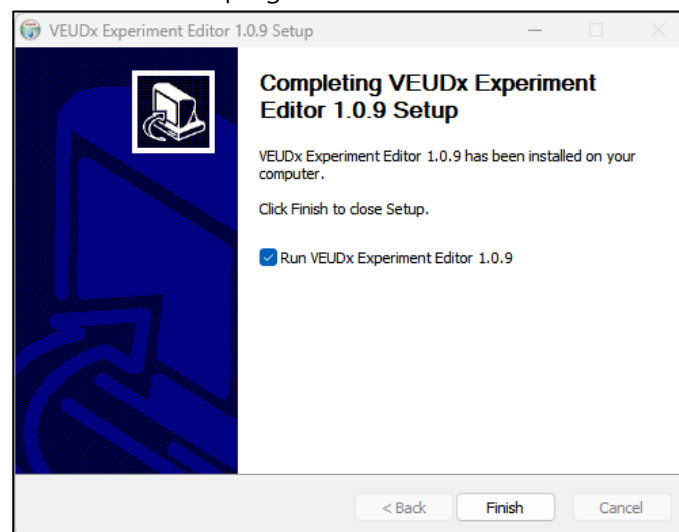

## 2. Term definition

### 2.1 ITEM

|                        |   |
|------------------------|---|
| <b>2.1.1 ITEM</b>      | ITEM consists of ITEM name( ex) TBI, Neurology ), Marker name, Pixel Cut, Experiment Protocol, etc. |
| <b>2.1.2 Protocol</b>  | Protocol is a collection of experimental procedures (Steps).  |
| <b>2.1.3 Step</b>      | This is the procedure for each well. (ex) Well 6 Washing 1 min )                                    |
| <b>2.1.3 ITEM file</b> | ITEM File is created as VEUDxITEM_ITEM_name.zip file.<br>(ex, VEUDxITEM_TBI.zip )                   |

### 2.2 LOT

|                       |   |
|-----------------------|---|
| <b>2.1.1 LOT</b>      | Depending on the produced LOT, it consists of Made Date, Serial, Expire Date, etc. An ITEM file is required to create a LOT.                                    |
| <b>2.1.2 LOT file</b> | A LOT XML file that stores LOT information and a Barcode PDF file are created.<br>Ex)<br>VEUDxLOT_TBI_EZTB22111601.xml<br>VEUDxLOT_TBI_EZTB22111601_BarCode.pdf |


### 2.3 QC Material LOT

|                                   |   |
|-----------------------------------|---|
| <b>2.3.1 QC Material LOT</b>      | Depending on the produced QC Material LOT, it consists of Made Date, Serial, Expire Date, etc. An ITEM file is required to create a QC Material LOT.                              |
| <b>2.3.2 QC Material LOT file</b> | QC Material LOT XML file that stores LOT information and a Barcode PDF file are created.<br>Ex)<br>VEUDx_QC_LOT_TBI-assay_230921001.xml<br>VEUDx_QC_LOT_TBI_230921001_BarCode.pdf |

## 3. Start

### 3.1 Start screen

**3.1.1 Start screen** On the start screen, you can select 'Make ITEM', 'Make LOT', 'Make QC Material LOT' or 'SW Manual'.

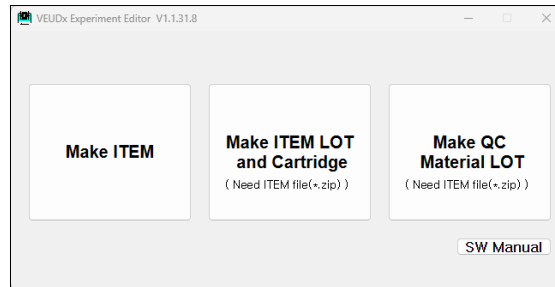

## 4. Make ITEM

### 4.1 ITEM creation start screen

#### 4.1.1 Start screen

ITEM Open, save and edit are possible.

- In the ITEM Summary tab, you can edit the item name, maker name by RSMP length, unit, CutOff, QC High/Low Range, QC Replication, pixel cut, experiment time, Dilution Factor, Fluorescence Exposure Time, Optical photography only, etc.

- You can edit steps in the Protocol tab

| Length | Marker | Use ExtMFI               | Count | Unit  | CutOff | QC Conc.<br>QC1 | QC Conc.<br>QC2 | LOD  | minMFI | QC1 min | QC1 Max | QC2 min | QC2 Max |
|--------|--------|--------------------------|-------|-------|--------|-----------------|-----------------|------|--------|---------|---------|---------|---------|
| 130    |        | <input type="checkbox"/> |       | pg/ml |        |                 |                 |      |        |         |         |         |         |
| 200    | UCH-L1 | <input type="checkbox"/> |       | pg/ml | 138    | 200             | 50              | 0.78 | 1500   | 5000    |         | 3000    |         |
| 250    |        | <input type="checkbox"/> |       | pg/ml |        |                 |                 |      |        |         |         |         |         |
| 300    | GFAP   | <input type="checkbox"/> |       | pg/ml | 24     | 100             | 25              | 0.39 | 1500   | 5000    |         | 3000    |         |
| 350    |        | <input type="checkbox"/> |       | pg/ml |        |                 |                 |      |        |         |         |         |         |
| 400    |        | <input type="checkbox"/> |       | pg/ml |        |                 |                 |      |        |         |         |         |         |
| 450    |        | <input type="checkbox"/> |       | pg/ml |        |                 |                 |      |        |         |         |         |         |
| 500    |        | <input type="checkbox"/> |       | pg/ml |        |                 |                 |      |        |         |         |         |         |

Total Estimated Time (Min.) 58  
 \* RT 6 samples assumptions  
 Reacting Time Only 49  
 Pixel Cut (%) Bottom 25 Top 5  
 Dilution Factor 4  
 Buffer Well None  
 QC Replication Not Specified  
 Fluorescence Exposure Time1(ms) 500  
 Fluorescence Exposure Time2(ms) 0  
☐ Use Extended MFI  
☐ Make RSMP QC Data  
☐ Optical photography only  
☐ Optical + Demagnetization only

Open ITEM Save ITEM



- You must input information about the name, unit, and RSMP Length of the markers to be used for the item.
- Based on the input information, the fluorescence image is detected by RSMP Length and the result is calculated.

<VEUDx Analysis Result Screen>

(Plus or minus 15% based on the center)

<VEUDx QC Result Screen>

## 4.2 Edit ITEM Summary

**4.2.1 Open ITEM** Press the “Open ITEM” button to select the previously created ITEM file.  
ITEM files are in \*.zip. (ex, VEUDxITEM\_TBI.zip )

VEUDx Experiment Editor V1.1.31.8 - Make ITEM

ITEM Summary Protocol Home Setting

ITEM Name TBI-assay Initialize

| Length | Marker | Use ExtMFI               | Count | Unit  | CutOff | QC Conc.<br>QC1 QC2 | LOD  | minMFI | QC1 min | QC1 Max | QC2 min | QC2 Max |
|--------|--------|--------------------------|-------|-------|--------|---------------------|------|--------|---------|---------|---------|---------|
| 130    |        | <input type="checkbox"/> |       | pg/ml |        |                     |      |        |         |         |         |         |
| 200    | UCH-L1 | <input type="checkbox"/> |       | pg/ml | 138    | 200 50              | 0.78 | 1500   | 5000    |         | 3000    |         |
| 250    |        | <input type="checkbox"/> |       | pg/ml |        |                     |      |        |         |         |         |         |
| 300    | GFAP   | <input type="checkbox"/> |       | pg/ml | 24     | 100 25              | 0.39 | 1500   | 5000    |         | 3000    |         |
| 350    |        | <input type="checkbox"/> |       | pg/ml |        |                     |      |        |         |         |         |         |
| 400    |        | <input type="checkbox"/> |       | pg/ml |        |                     |      |        |         |         |         |         |
| 450    |        | <input type="checkbox"/> |       | pg/ml |        |                     |      |        |         |         |         |         |
| 500    |        | <input type="checkbox"/> |       | pg/ml |        |                     |      |        |         |         |         |         |

Total Estimated Time (Min.) 58 QC Replication Not Specified

\* RT 6 samples assumptions

Reacting Time Only 49

Fluorescence Exposure Time1(ms) 500

Fluorescence Exposure Time2(ms) 0

Pixel Cut (%) Dilution Factor 4

Bottom 25

Top 5 Buffer Well None

Use Extended MFI

Make RSMP QC Data

Optical photography only

Optical + Demagnetization only

Open ITEM Save ITEM

**4.2.2 Edit ITEM Summary** Editable below items

- Item name
- Marker name used by RSMP length
- Unit used by RSMP length
- CutOff for Result(Positive/Negative)
- QC High/Low Range



If you want to use a different concentration unit, you can add a concentration unit in 'Setting'

- Pixel Cut Bottom, Top (%)  
Remove noise caused by saturation when obtaining MFI.  
Defaults to Bottom 25% Top 5%.
- Dilution Factor  
The calculated concentration value is multiplied by the 'dilution factor' value.



- Fluorescence Exposure Time

The default value of Fluorescence Exposure is 500ms.(1~1000)

If the MFI value of the item you are using is high or low, adjust the 'Fluorescence Exposure' value.

(High MFI values can affect fluorescence saturation)

- Optical photography only

Check if you only want to image capture and Detecting (RSMP is demagnetized)

- Optical + Demagnetization only


Check if you only want to image capture and Detecting (RSMP is not demagnetized)

(RSMP should be in imaging well 8)

- Make RSMP QC Data

If you select the 'Make RSMP QC Data' Option, you can acquire MFI, CV(%) values for each RSMP and for the entire RSMP.

(Used to check RSMP Ab coupling QC status)

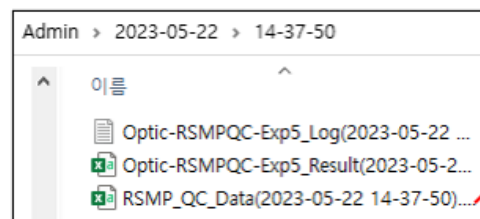

<Select Make RSMP QC Data Option>


| Cartridge1 |         |         |         |
|------------|---------|---------|---------|
| RSMP       | MFI     | STDEV   | CV(%)   |
| 1          | 27793.9 | 7670.2  | 27.5967 |
| 2          | 25359.5 | 6268.25 | 24.7176 |
| 3          | 32523.1 | 9540.35 | 29.3341 |
| 4          | 29816.1 | 8989.08 | 30.1484 |
| 5          | 26123.3 | 6210.1  | 23.7722 |
| 6          | 30254.4 | 7411.64 | 24.4977 |
| 7          | 25311.5 | 8759.78 | 34.6079 |
| 8          | 22506.7 | 7082.12 | 31.4668 |
| 9          | 27966.5 | 7981.9  | 28.541  |
| 10         | 29312   | 6347.73 | 21.6557 |
| 11         | 29546.5 | 9441.77 | 31.9556 |
| 12         | 32746.4 | 15021.4 | 45.872  |
| 13         | 29371   | 12198.2 | 41.5313 |
| 14         | 25405.5 | 11669.9 | 45.9346 |
| 15         | 28429   | 11555.8 | 40.648  |
| 16         | 20178.1 | 4986.54 | 24.7126 |
| 17         | 24910.8 | 5694.06 | 22.8578 |
| 18         | 24262   | 8428.5  | 34.7395 |
| 19         | 31228.1 | 7274.75 | 23.2955 |
| 20         | 29123.8 | 5586.85 | 19.1831 |
| 21         | 26226.2 | 6798.9  | 25.9241 |
| 22         | 28735.2 | 6419.81 | 22.3413 |
| 23         | 26659.1 | 8307.48 | 31.1618 |
| 24         | 29073.8 | 6705.22 | 23.0628 |
| 25         | 29002.5 | 5813.24 | 20.044  |
| 26         | 33350.8 | 10017.8 | 30.0377 |
| 27         | 31911.5 | 9865.93 | 30.9165 |
| 28         | 32402.5 | 9230.65 | 28.4875 |
| 29         | 28249.9 | 7301.46 | 25.846  |
| 30         | 29633.6 | 6110.93 | 20.6216 |
| Total      | MFI     | STDEV   | CV(%)   |
|            | 28247.1 | 3057.53 | 10.8242 |

<RSMP QC Data Result(.csv) file>

### 4.3 Edit Protocol


#### 4.3.1 Protocol tab

If you select the Protocol tab, you can edit the Step.

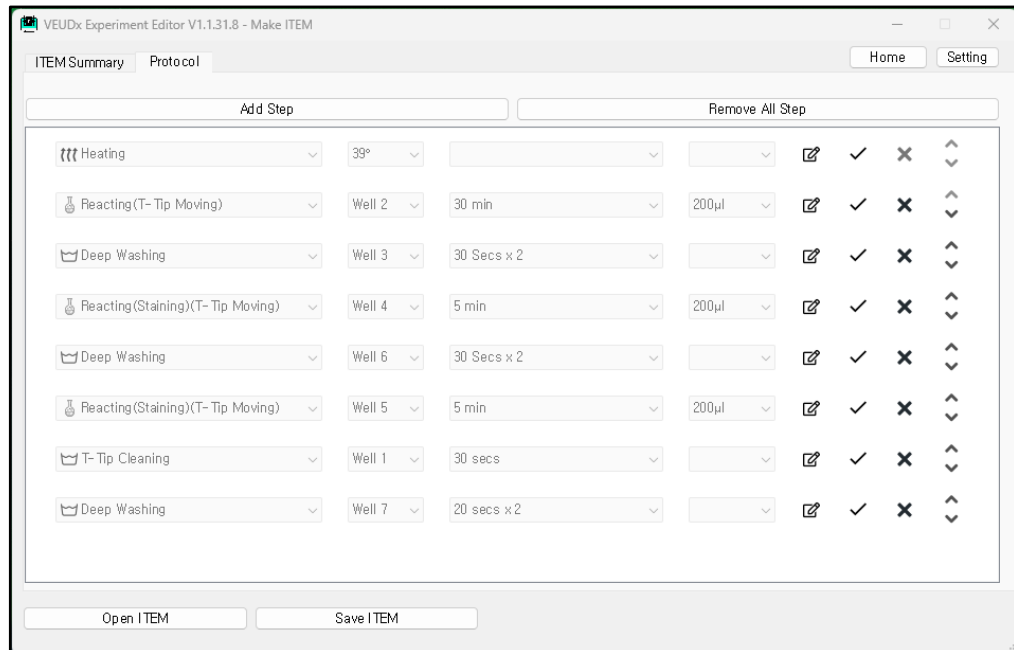

#### 4.3.2 Add Step

Click the "Add Step" box at the top and add an item.

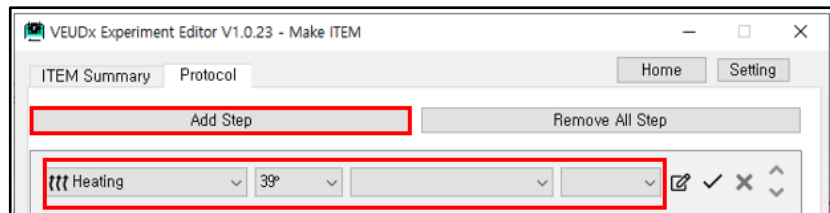

- Selection by Category

.Heating : 30° / 35° / 37° / 38° / 39° / X / RT selectable

.Reacting(Staining)(T-Tip Moving)

: Well2~7 / Time(3 min, 5 min) / T-Tip depth (150µl, 200µl) selectable

\*Shuffle every 2min

.Reacting(T-Tip Moving)

: Well2~7 / Time (1~480 min) / T-Tip depth (150µl, 200µl) selectable

\*Shuffle every 5min

.Reacting Shuffle(T-Tip Moving)

: Well2~7 / Shuffle Time (1,2,5 min) / Time( 1~60 )

.Reacting(Staining)(M-Bar Moving) (Not Recommend)

: Well2~7 / Time (3 min, 5 min) / T-Tip depth (150μl, 200μl) selectable

.Reacting(M-Bar Moving) (Not Recommend)


: Well2~7 / Time (1~480 min) / T-Tip depth (150μl, 200μl) selectable

.Washing : Well1~7 / Time ( 10, 20, 30 secs, 1 min, 30 secs x4 ) selectable

.Deep Washing : Well1~7 / Time ( 20 secs, 30 secs, 20 secs x 2, 1 min,  
30 secsX2, 1 min+30 secs X 3) selectable

.T-Tip Cleaning : Well1~7 / Time ( 20, 30, 45, 90 secs) selectable

'T-tip Cleaning' use after Staining Reaction to prevent PE buffer from remaining on the T-Tip into the Imaging Well.

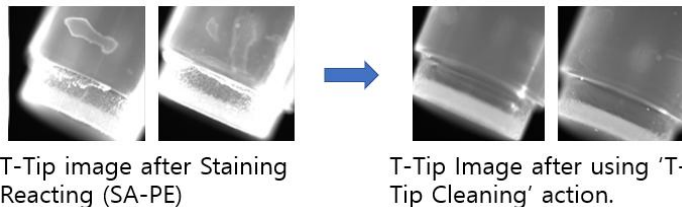

Heating can be selected only once and must be located only in the first step.

#### 4.3.3 Edit Step



- : Change edit Step mode
- : Save Step
- : Delete Step
- : Change Step Order


## 4.4 Save ITEM and Installation

### 4.4.1 Save ITEM to PC

Click the "Save ITEM" button to create an ITEM zip file.

\*For how to save only protocol (script) for development reference, refer to the setting section.

- 4.4.2 ITEM installation on VEUDx equipment**
1. Copy the ITEM file created above to an external USB memory
  2. Run VEUDx equipment
  3. Admin Login (Initial Admin Password: 0000)
  4. Click Settings
  5. Click Information
  6. Click Item
  7. Mounting on an external USB memory device
  8. Select Update

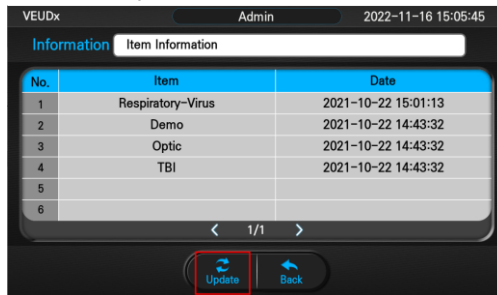


When you press "Update" button, you can check the updateable ITEM list.  
(Update file must be placed in the USB Root folder)



9. Select ITEM to update
10. ITEM update complete
11. Restart after shutting down the equipment



When you press each item, you can check detailed information such as each marker name and unit.

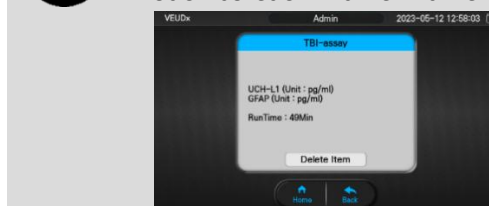


## 5. Make LOT

### 5.1 Select ITEM file


#### 5.1.1 Select ITEM file

Select ITEM file to make LOT.

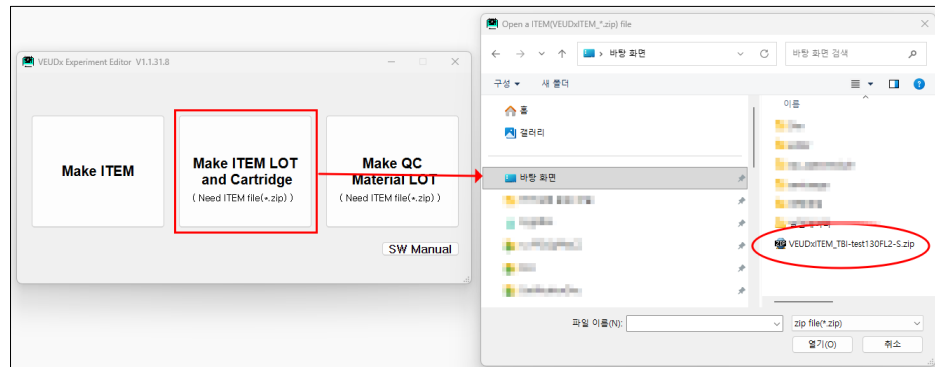

### 5.2 Make LOT

#### 5.2.1 Start

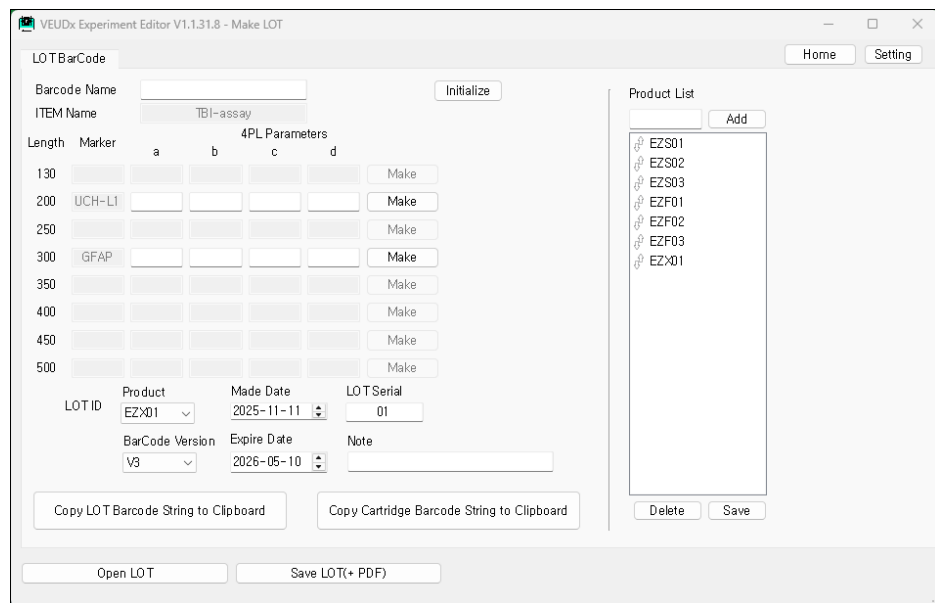

The ITEM name and Marker name read from the ITEM file are displayed.

#### 5.2.2 Enter LOT information

- 4PL Parameters (a, b, c, d)
  - LOT creation date
  - LOT Serial
  - LOT Expire Date can be entered.
- \*Barcode Name, Note can be used as a reference.


#### 5.2.3 Calculate 4PL Parameters

- Select "Make" button for each Maker.

| Length | Marker | 4PL Parameters |   |   |   |
|--------|--------|----------------|---|---|---|
|        |        | a              | b | c | d |
| 200    | 200um  |                |   |   |   |

**Make**

- 4PL calculation screen

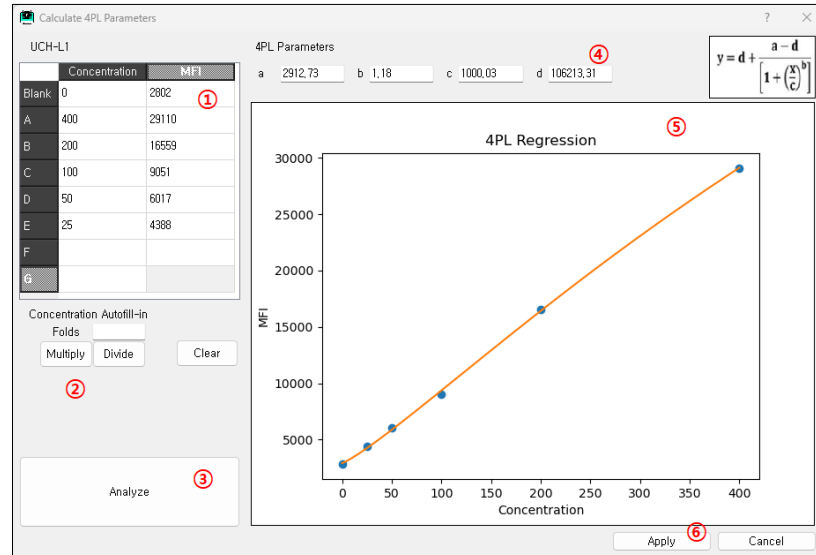

- ① Concentration, MFI input  
(Copy/Paste available in Excel/Notepad)
- ② Automatically multiply and divide as much as Folds(multiples) based on A concentration
- ③ Calculate 4PL parameters and graphs according to the input concentration and MFI values
- ④ Calculated graph
- ⑤ Apply to the LOT edit screen by selecting "Apply" button

VEUDx Experiment Editor V1.0.23 - Make LOT

Home Setting

Barcode Name

ITEM Name QC-ex5-W3-H37-2

Initialize

4PL Parameters

| Length | Marker | a       | b    | c      | d        |
|--------|--------|---------|------|--------|----------|
| 200    | 200um  | 2856.50 | 1.20 | 885.27 | 97090.44 |
| 250    |        |         |      |        |          |
| 300    | 300um  |         |      |        |          |
| 350    |        |         |      |        |          |
| 400    | 400um  |         |      |        |          |
| 450    |        |         |      |        |          |
| 500    | 500um  |         |      |        |          |

Make

LOT ID

Maker EZ ITEM Abbr. QC Made date 2023-05-16 LOT Serial 01

Bar Code Version V1 Expire Date 2023-11-12

Note

Open LOT Save LOT(+ PDF)

### 5.2.4 Copy Barcode string to Clipboard

When you press the "Copy LOT Barcode string to Clipboard" Button, the LOT Barcode string below will be copied to the clipboard.

Ex)"VEUDx-L/V3/230927001/TBI-

V2/240325/1\_a2\_b3\_c4\_d5\_x65535\_n0/3\_a6\_b7\_c8\_d9\_x65535\_n0"

When you press the "Copy Cartridge Barcode string to Clipboard" Button, the Cartridge Barcode string below is copied to the clipboard.

ex)"VEUDx-C/V3/230927001"

### 5.2.3 Save LOT

Click the "Save LOT(+PDF)" button


Two LOT XML files and Barcode PDF files are created.

- LOT XML file (ex, VEUDxLOT\_TBI\_EZTB22111601.xml) Used to save work

When using VEUDx equipment, it can be used when there is no Barcode Reader equipment.

- Barcode PDF file (ex, VEUDxLOT\_TBI\_EZTB22111601\_BarCode.pdf)

It is used for reading VEUDx equipment LOT/Script Barcode by printing it out on paper.

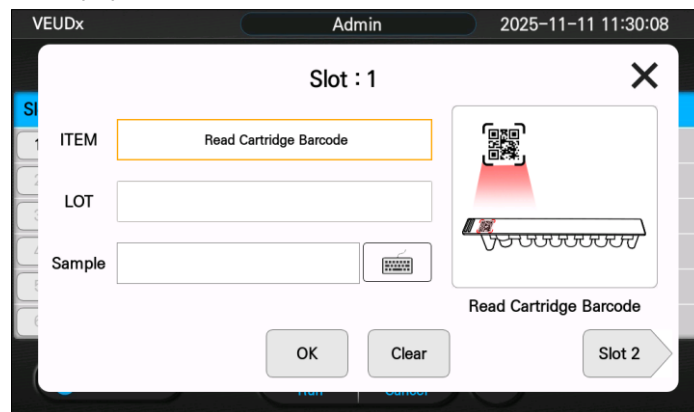

< VEUDx equipment LOT reading screen >



< PDF file for LOT Barcode output >

**5.2.4 Open LOT** Editing is possible by selecting the saved LOT XML by pressing the "Open LOT" button.




## 6. Make QC Material LOT

### 6.1 Select ITEM file


#### 6.1.1 Select ITEM file

Select ITEM file to make QC Material LOT.

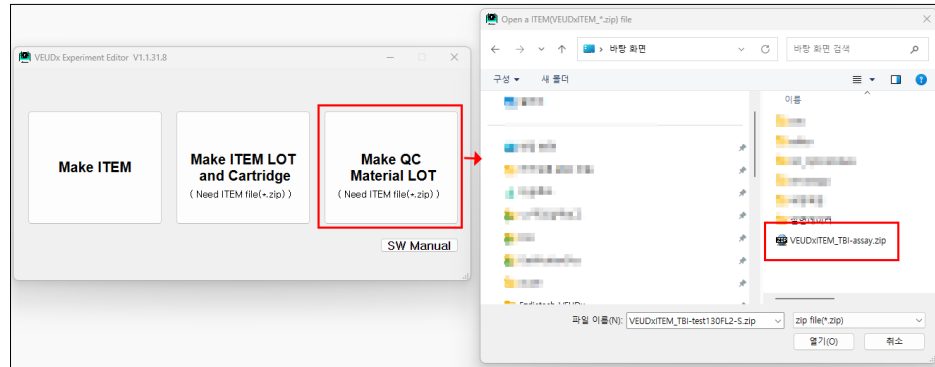

### 6.2 Make QC Material LOT

#### 6.2.1 Start

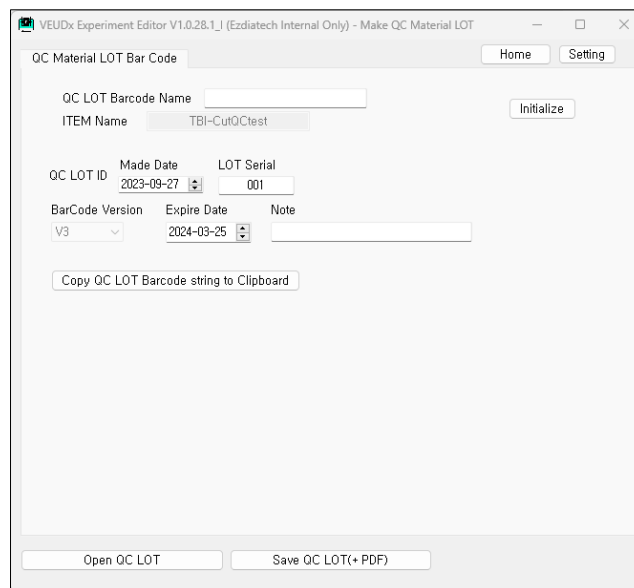

The ITEM name read from the ITEM file are displayed.

#### 6.2.2 Enter QC Material LOT information

- QC Material LOT creation date
- QC Material LOT Serial
- QC Material LOT Expire Date can be entered.
- \*Barcode Name, Note can be used as a reference.

#### 6.2.3 Copy Barcode string

When you press the "Copy QC LOT Barcode string to Clipboard" Button, the QC LOT Barcode string below will be copied to the clipboard.

**to Clipboard** Ex)" VEUDx-Q/V3/230927001/TBI-CutQCtest/240325"


**6.2.4 Save QC** Click the "Save QC LOT(+PDF)" button

**Material LOT** Two QC LOT XML files and Barcode PDF files are created.

- QC Material LOT XML file (ex, VEUDx\_QC\_LOT\_TBI-assay\_230922001.xml)  
used to save work

- QC Material Barcode PDF file  
(ex, VEUDx\_QC\_LOT\_TBI-assay\_230922001\_BarCode.pdf)

It is used for reading VEUDx equipment QC Material Barcode by printing it out on paper.

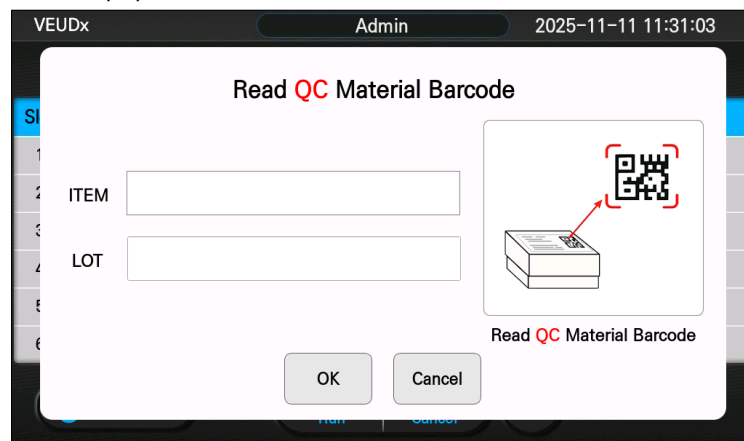

< VEUDx equipment QC Material LOT reading screen >



< PDF file for QC Material LOT Barcode output >


**6.2.5 Open QC** Editing is possible by selecting the saved LOT XML by pressing the "Open  
**Material LOT** QC LOT" button.

## 6. Settings

### 6.1 Open Settings

#### 6.1.1 Open Settings


'Press the 'Settings' button.

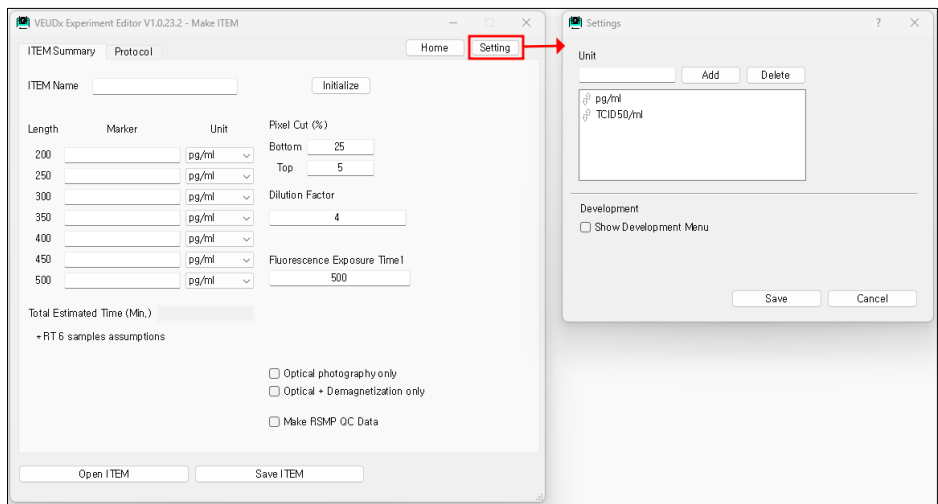

### 6.2 Unit

#### 6.2.1 Unit


Units to be displayed in the Unit Combo List of the ITEM Summary tab can be edited.

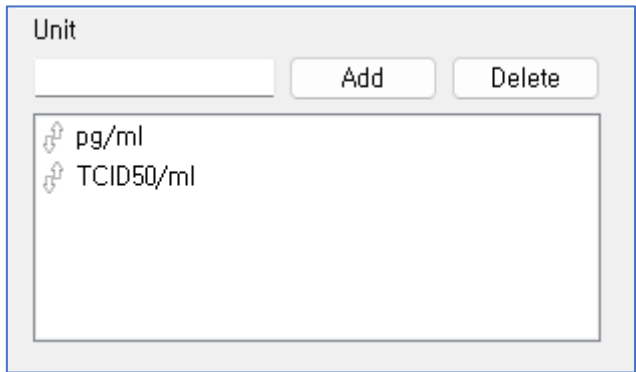

### 6.3 Development menu

#### 6.3.1 Development menu

If you select the development menu, you can open/save only the protocol and display comments for analysis.

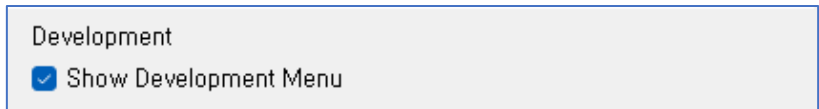

#### 6.4.2 Save protocol file

If you press the 'Save Protocol' button, only the Protocol (Script) file can be saved separately.

This is a development file and **cannot** be installed on a machine.

When "Add Comment", comments are displayed on the Protocol (Script).

< Protocol display according to annotation options >