
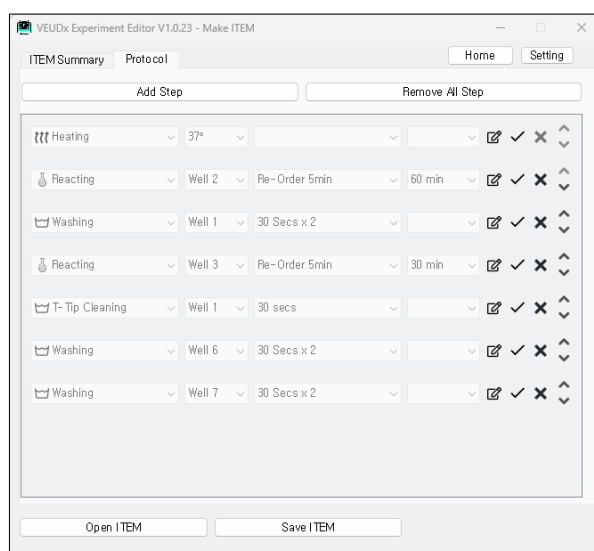


Operation Manual

VEUDx Experiment Editor

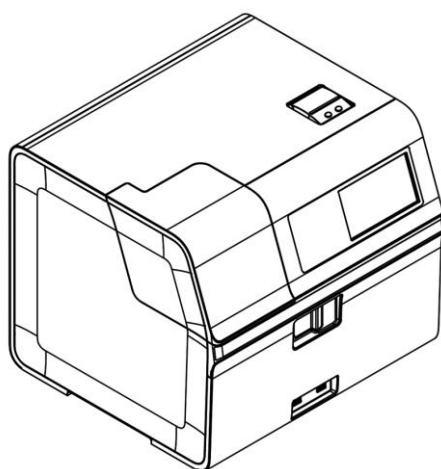

UM-VEUDx-1.0

Index

1. Installation.....	3
2. Term definition.....	5
2.1 ITEM	5
2.2 LOT	5
3. Start.....	6
3.1 Start screen.....	6
4. Make ITEM	6
4.1 ITEM creation start screen	6
4.2 Edit ITEM Summary	8
4.3 Edit Protocol	11
4.4 Save ITEM and Installation	12
5. Make LOT	14
5.1 Select ITEM file	14
5.2 Make LOT	14
6. Settings.....	17
6.1 Open Settings.....	17
6.2 Unit.....	18
6.3 Development menu.....	18

1. Installation


1.1 Installation

1.1.1 Execution

Execute VEUDx Experiment Editor Setup V1.x.x.exe on PC.


1.1.2 Language selection

Choose the language you want to use

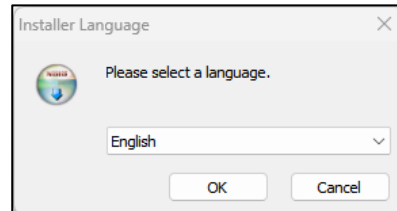


1.1.3 Start installation

Check the contents and click 'Next'.

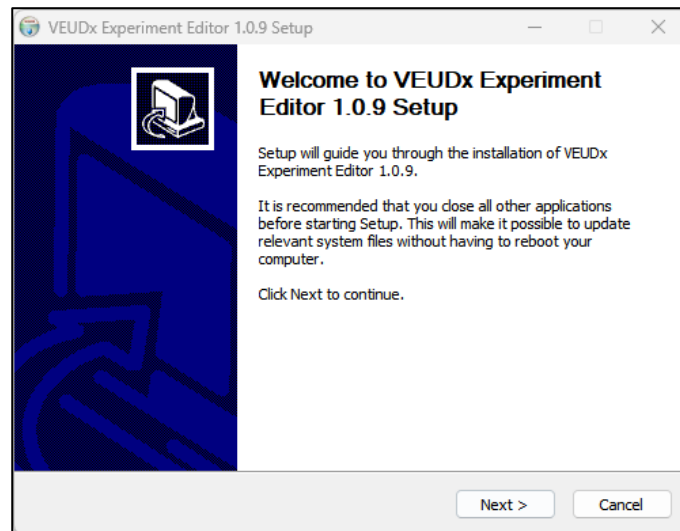


1.1.4 Check License

Read the license agreement and click 'I Agree'.

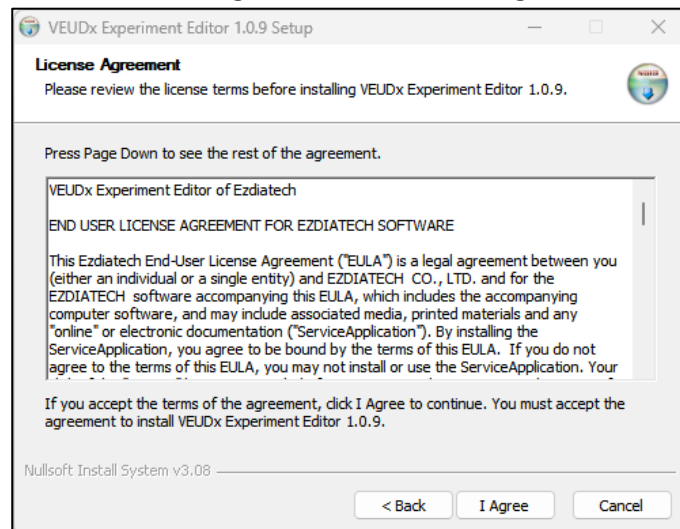


1.1.5 Installation location

Check the installation path and click 'Install'.

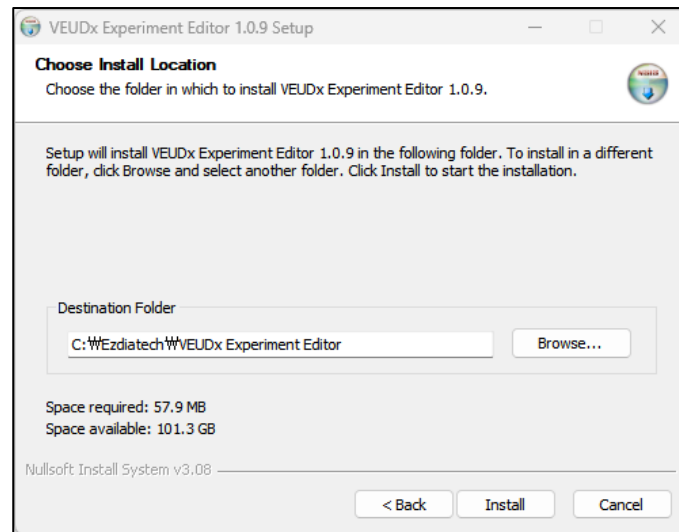

1.1.6 Installation completion

Check to run the program and click 'Finish'.

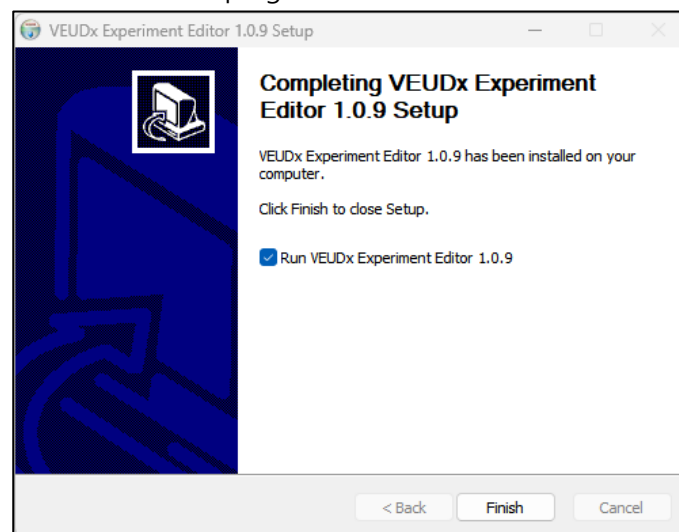

2. Term definition

2.1 ITEM

- 2.1.1 ITEM** ITEM consists of ITEM name(ex) TBI, Neurology), Marker name, Pixel Cut, Experiment Protocol, etc.
- 2.1.2 Protocol** Protocol is a collection of experimental procedures (Steps).
- 2.1.3 Step** This is the procedure for each well. (ex) Well 6 Washing 1 min)
- 2.1.3 ITEM file** ITEM File is created as VEUDxITEM_ITEM_name.zip file.
(ex, VEUDxITEM_TBI.zip)


2.2 LOT

- 2.1.1 LOT** Depending on the produced LOT, it consists of calibration, concentration cut-off value, etc. An ITEM file is required to create a LOT.
- 2.1.2 LOT file** A LOT XML file that stores LOT information and a Barcode PDF file are created.
Ex)
VEUDxLOT_TBI_EZTB22111601.xml
VEUDxLOT_TBI_EZTB22111601_BarCode.pdf

3. Start

3.1 Start screen

3.1.1 Start screen On the start screen, you can select 'Make ITEM', 'Make LOT', or 'SW Manual'.

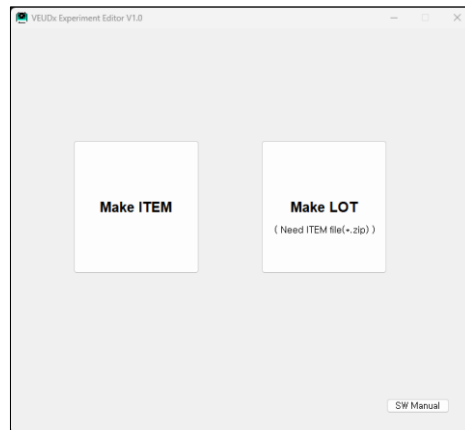

4. Make ITEM

4.1 ITEM creation start screen

4.1.1 Start screen ITEM Open, save and edit are possible.

- In the ITEM Summary tab, you can edit the item name, maker name by RSMP length, unit, pixel cut, experiment time, Dilution Factor, Fluorescence Exposure Time, Optical photography only, etc.
- You can edit steps in the Protocol tab

Length	Marker	Unit	Pixel Cut (%)
200		pg/ml	Bottom 25
250		pg/ml	Top 5
300		pg/ml	Dilution Factor 4
350		pg/ml	
400		pg/ml	
450		pg/ml	Fluorescence Exposure Time 1 500
500		pg/ml	



- You must input information about the name, unit, and RSMP Length of the markers to be used for the item.
- Based on the input information, the fluorescence image is detected by RSMP Length and the result is calculated.

Length	Marker	Unit
200	UCH-L1	pg/ml
250		TCID50/ml
300	GFAP	pg/ml
350		TCID50/ml
400		TCID50/ml
450		TCID50/ml
500		TCID50/ml

TBI-assay	Slot1	Slot2	Slot3	Slot4	Slot5	Slot6
Sample ID : EZ-R-001	UCH-L1 : 399.722(pg/ml) GFAP : 199.948(pg/ml)					

<VEUDx Result Screen>

ResearcherID : Admin
 AnalysisDate : 2023-05-12 14:50:27
 RunningTime : 00:48:37
 ITEM : TBI-assay
 LOT ID : EZTB23051201

MFI value are converted to Concentrations via the 4PL Formula

MFI value for each cartridge and marker

#Sample	Conc.		DF		MFI		Well Bg Brightness	Well Bg CV(%)
	UCH-L1	GFAP	UCH-L1	GFAP	UCH-L1	GFAP		
EZ-R-001	399.722	199.948	1598.89	799.792	21638.7	37300.6	3697.77	41.5642
EZ-R-002	201.038	100.511	804.152	402.044	12333.3	16572.3	3463.81	26.0463
EZ-R-003	91.8157	45.9618	367.263	183.847	5409.47	7167.85	3344.48	21.123
EZ-R-004	74.8248	31.1242	299.299	124.497	4744.2	5598.43	3418.72	15.8736
EZ-R-005	51.1242	23.129	204.497	92.510	4107.44	4987.54	3457.1	21.2582

pg/ml pg/ml

ITEM(LOT ID)
 TBI-assay(EZTB23051201)

Marker	Cal_a	Cal_b	Cal_c	Cal_d
UCH-L1	3706.14	2.57	245.62	26768.45
GFAP	4225.27	2.01	226.57	79823.25

Standard Curve data (4PL parameters) by Marker

Pixel Top cut : 5
 Pixel Bottom cut : 25
 FL_LED : 361

<VEUDx Result(.csv) file>

4.2 Edit ITEM Summary

4.2.1 Open ITEM

Press the "Open ITEM" button to select the previously created ITEM file.
 ITEM files are in *.zip. (ex, VEUDxITEM_TBI.zip)

VEUDx Experiment Editor V1.0.23.2 - Make ITEM

ITEM Summary Protocol Home Setting

ITEM Name Initialize

Length	Marker	Unit	Pixel Cut (%)
200	UCH-L1	pg/ml	Bottom 25
250		pg/ml	Top 5
300	GFAP	pg/ml	Dilution Factor 4
350		pg/ml	Fluorescence Exposure Time1 500
400		pg/ml	
450		pg/ml	
500		pg/ml	

Total Estimated Time (Min.)

* RT 6 samples assumptions

☐ Optical photography only
☐ Optical + Demagnetization only
☐ Make RSMP QC Data

Open ITEM Save ITEM

4.2.2 Edit Editable below items**ITEM**

- Item name

Summary

- Marker name used by RSMP length
- Unit used by RSMP length
- Total experiment time
- Unit



If you want to use a different concentration unit, you can add a concentration unit in 'Setting'

- Pixel Cut Bottom, Top (%)

Remove noise caused by saturation when obtaining MFI.

Defaults to Bottom 25% Top 5%.

- Dilution Factor

The calculated concentration value and the 'Dilution Factor' value are multiplied and displayed separately in a VEUDx Result(.csv) File.

ResearcherID : Admin								
AnalysisDate : 2023-05-12 14:50:27								
RunningTime : 00:48:37								
ITEM : TBI-assay								
LOT ID : EZTB23051201								
#.Sample	Conc.		*DF		MFI		Well Bg Brightness	Well Bg CV(%)
	UCH-L1	GFAP	UCH-L1	GFAP	UCH-L1	GFAP		
EZ-R-001	399.722	199.948	1598.89	799.792	21638.7	37300.6	3697.77	41.5642
EZ-R-002	201.038	100.511	804.152	402.044	12333.3	16572.3	3463.81	26.0463
EZ-R-003	91.8157	45.9618	367.263	183.847	5409.47	7167.85	3344.48	21.123
EZ-R-004	74.8248	31.1242	299.299	124.497	4744.2	5598.43	3418.72	15.8736
EZ-R-005	51.1242	23.129	204.497	92.516	4107.44	4987.54	3457.1	21.2582
	pg/ml	pg/ml						
ITEM(LOT ID)	Marker	Cal_a	Cal_b	Cal_c	Cal_d			
TBI-assay(EZTB23051201)	UCH-L1	3706.14	2.57	245.62	26768.45			
	GFAP	4225.27	2.01	226.57	79823.25			
Pixel Top cut : 5								
Pixel Bottom cut : 25								
FL_LED : 361								

<VEUDx Result(.csv) file>

- Fluorescence Exposure Time

The default value of Fluorescence Exposure is 500ms.(1~1000)

If the MFI value of the item you are using is high or low, adjust the 'Fluorescence Exposure' value.

(High MFI values can affect fluorescence saturation)

- Optical photography only

Check if you only want to image capture and Detecting

(RSMP is demagnetized)

- Optical + Demagnetization only

Check if you only want to image capture and Detecting


(RSMP is not demagnetized)

(RSMP should be in imaging well 8)

- Make RSMP QC Data

If you select the 'Make RSMP QC Data' Option, you can acquire MFI, CV(%) values for each RSMP and for the entire RSMP.

(Used to check RSMP Ab coupling QC status)

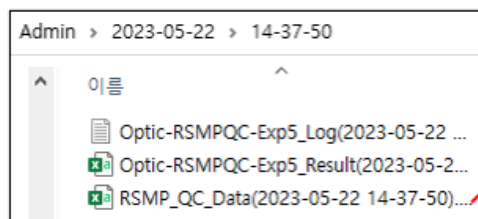

<Select Make RSMP QC Data Option>

Cartridge1			
RSMP	MFI	STDEV	CV(%)
1	27793.9	7670.2	27.5967
2	25359.5	6268.25	24.7176
3	32523.1	9540.35	29.3341
4	29816.1	8989.08	30.1484
5	26123.3	6210.1	23.7722
6	30254.4	7411.64	24.4977
7	25311.5	8759.78	34.6079
8	22506.7	7082.12	31.4668
9	27966.5	7981.9	28.541
10	29312	6347.73	21.6557
11	29546.5	9441.77	31.9556
12	32746.4	15021.4	45.872
13	29371	12198.2	41.5313
14	25405.5	11669.9	45.9346
15	28429	11555.8	40.648
16	20178.1	4986.54	24.7126
17	24910.8	5694.06	22.8578
18	24262	8428.5	34.7395
19	31228.1	7274.75	23.2955
20	29123.8	5586.85	19.1831
21	26226.2	6798.9	25.9241
22	28735.2	6419.81	22.3413
23	26659.1	8307.48	31.1618
24	29073.8	6705.22	23.0628
25	29002.5	5813.24	20.044
26	33350.8	10017.8	30.0377
27	31911.5	9865.93	30.9165
28	32402.5	9230.65	28.4875
29	28249.9	7301.46	25.846
30	29633.6	6110.93	20.6216
Total	MFI	STDEV	CV(%)
	28247.1	3057.53	10.8242


<RSMP QC Data Result(.csv) file>

4.3 Edit Protocol

4.3.1


If you select the Protocol tab, you can edit the Step.

Protocol tab

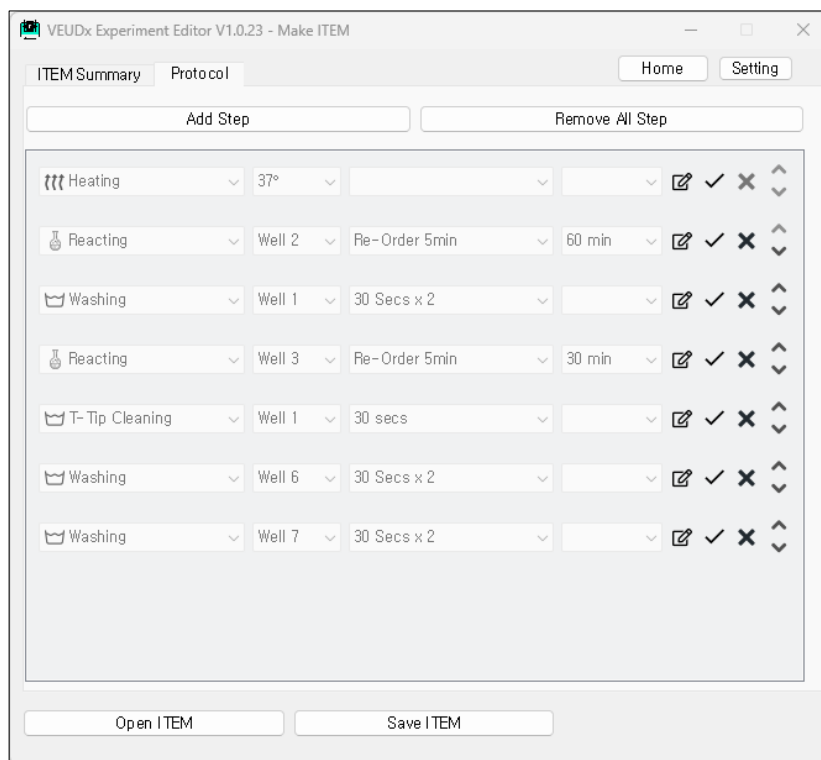

4.3.2 Add Step

Click the "Add Step" box at the top and add an item.

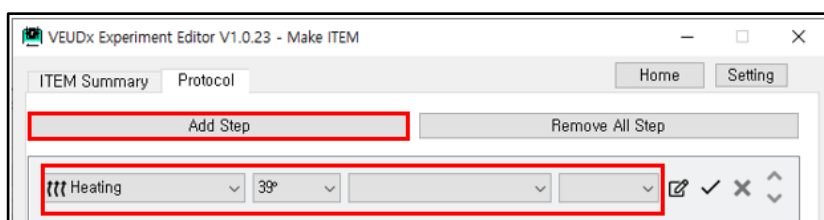

- Selection by Category

.Heating : 30° / 35° / 37° / 38° / 39° / X / RT selectable

.Reacting : Reaction Well (Well 2~5) selectable

- Re-Order 2min (Change the location of RSMP every 2 minutes.)

-> 1~30 min selectable


- Re-Order 5min (Change the location of RSMP every 5 minutes.)

-> 1~480 min selectable

.Washing: Well1~7 and time (30 secs, 20 secsX2, 30 secsX2) selectable


.T-Tip Cleaning: Well1~7 and time (30 secs) selectable



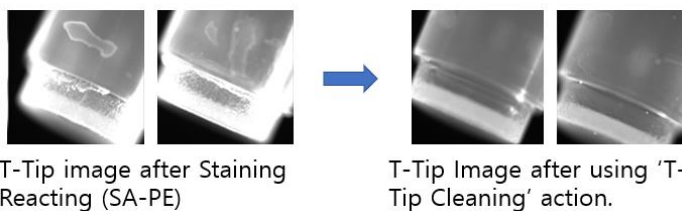

'T-tip Cleaning' use after Staining Reaction to prevent PE buffer from remaining on the T-Tip into the Imaging Well.

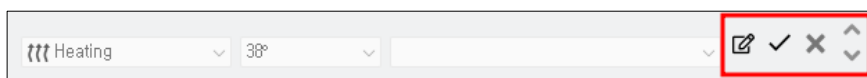

Heating can be selected only once and must be located only in the first step.

4.3.3 Edit Step



-  : Change edit Step mode
-  : Save Step
-  : Delete Step
-  : Change Step Order

4.4 Save ITEM and Installation

4.4.1 Save ITEM to PC


Click the "Save ITEM" button to create an ITEM zip file.

*For how to save only protocol (script) for development reference, refer to the setting section.

4.4.2 ITEM installation on VEUDx equipment

1. Copy the ITEM file created above to an external USB memory
2. Run VEUDx equipment
3. Admin Login (Initial Admin Password: 0000)
4. Click Settings
5. Click Information
6. Click Item
7. Mounting on an external USB memory device

8. Select Update

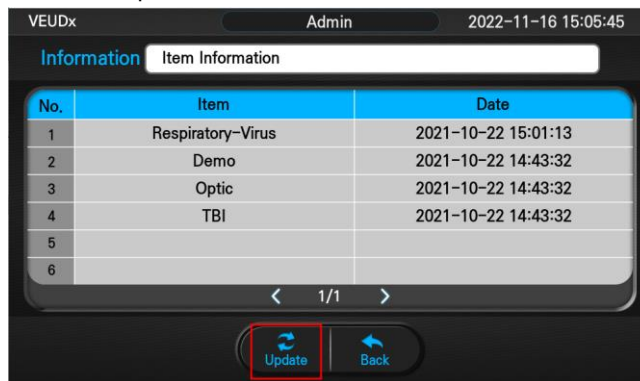

When you press "Update" button, you can check the updateable ITEM list.
(Update file must be placed in the USB Root folder)




9. Select ITEM to update

10. ITEM update complete

11. Restart after shutting down the equipment



When you press each item, you can check detailed information such as each marker name and unit.

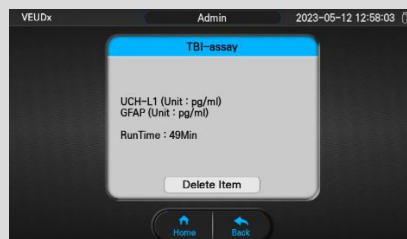


5. Make LOT

5.1 Select ITEM file

5.1.1 Select ITEM file

Select ITEM file to make LOT.

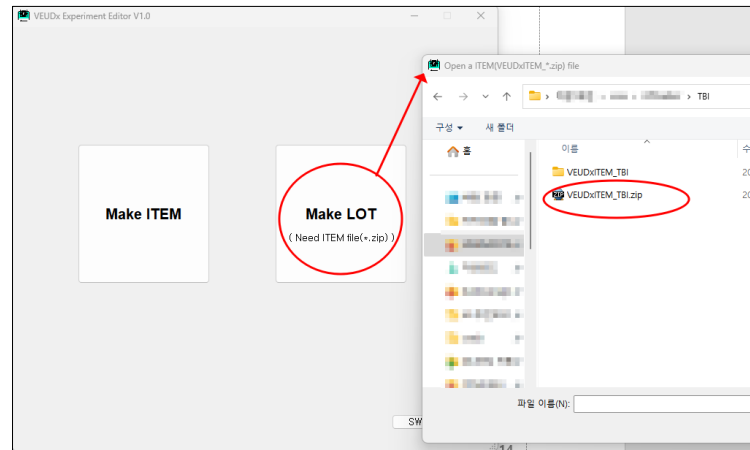

5.2 Make LOT

5.2.1 Start

The screenshot shows the 'VEUDx Experiment Editor V1.0.23 - Make LOT' window. It contains the following fields and controls:

- LOT Bar Code**: Home, Setting buttons
- Barcode Name**: Text input field
- ITEM Name**: Text input field with value 'QC-ex5-W3-H37-2', Initialize button
- 4PL Parameters**:

Length	Marker	a	b	c	d	
200	200um					Make
250						Make
300	300um					Make
350						Make
400	400um					Make
450						Make
500	500um					Make
- LOT ID**: Maker (EZ), ITEM Abbr. (QC), Made date (2023-05-16), LOT Serial (01)
- Bar Code Version**: VI, **Expire Date**: 2023-11-12
- Note**: Text input field
- Buttons**: Open LOT, Save LOT(+ PDF)

The ITEM name and Marker name read from the ITEM file are displayed.

5.2.2 Enter LOT information

- 4PL Parameters (a, b, c, d)
- LOT creation date
- LOT Serial
- LOT Expire Date can be entered.
- *Barcode Name, Note can be used as a reference.


5.2.3 Calculate 4PL Parameters

- Select "Make" button for each Maker.

Length	Marker	4PL Parameters			
		a	b	c	d
200	200um				

Make

- 4PL calculation screen

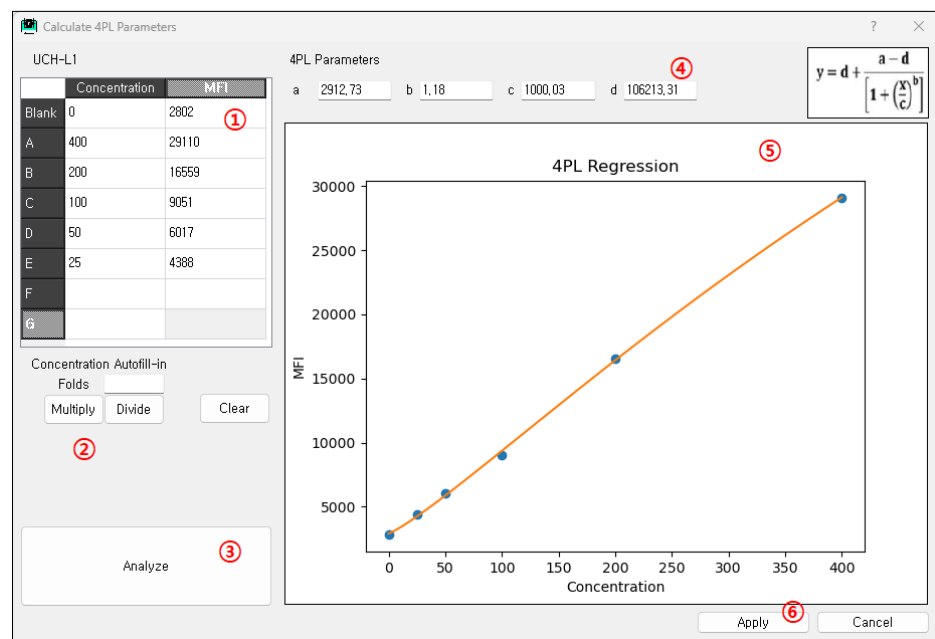

- ① Concentration, MFI input
(Copy/Paste available in Excel/Notepad)
- ② Automatically multiply and divide as much as Folds(multiples) based on A concentration
- ③ Calculate 4PL parameters and graphs according to the input concentration and MFI values
- ④ Calculated graph
- ⑤ Apply to the LOT edit screen by selecting "Apply" button

LOT Bar Code

Barcode Name:

ITCM Name:

4PL Parameters

Length	Marker	a	b	c	d	Make
200	200um	285.58	1.28	885.27	87080.44	Make
250						Make
300	300um					Make
350						Make
400	400um					Make
450						Make
500	500um					Make

LOT ID:

Bar Code Version:

Expire Date:

Buttons: Open LOT, Save LOT(+PDF)

5.2.3 Save LOT

Click the "Save LOT(+PDF)" button


Two LOT XML files and Barcode PDF files are created.

- LOT XML file (ex, VEUDxLOT_TBI_EZTB22111601.xml) Used to save work

When using VEUDx equipment, it can be used when there is no Barcode Reader equipment.

- Barcode PDF file (ex, VEUDxLOT_TBI_EZTB22111601_BarCode.pdf)

It is used for reading VEUDx equipment LOT/Script Barcode by printing it out on paper.

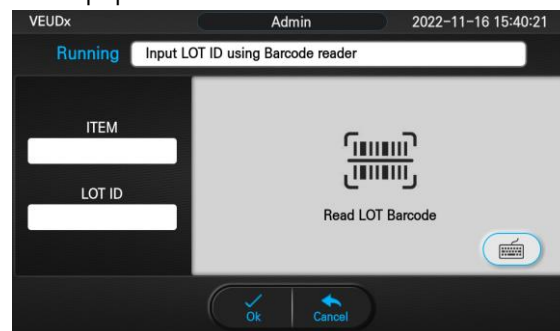

< VEUDx equipment LOT reading screen >

< PDF file for LOT Barcode output >

Editing is possible by selecting the saved LOT XML by pressing the "Open LOT" button.

'Press the 'Settings' button.



6.2 Unit

6.2.1 Unit

Units to be displayed in the Unit Combo List of the ITEM Summary tab can be edited.

6.3 Development menu

6.3.1

Development menu


If you select the development menu, you can open/save only the protocol and display comments for analysis.

6.4.2 Save protocol file

If you press the 'Save Protocol' button, only the Protocol (Script) file can be saved separately.

This is a development file and **cannot** be installed on a machine.

When "Add Comment", comments are displayed on the Protocol (Script).

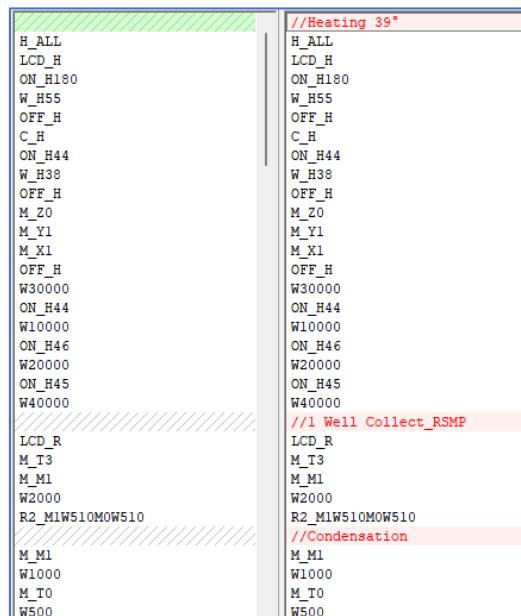

< Protocol display according to annotation options >